

# BladexT

The complete solution for an active lifestyle.






## **Blade**XT

BladeXT is the complete solution for active users. It is suitable for everyday use as well as running, cycling and other sporting activities, too. It has been designed so that wearers don't have to stop and change their prostheses between activities throughout the day.

BladeXT is lightweight, flexible and responsive, yet durable. It delivers outstanding performance, ground compliance, comfort and postural support from the twin toe spring and traction heel, while its C-shaped toe spring is primed for optimal energy response. The result is a sports blade that offers the user versatility, security and the confidence to enjoy an active life.

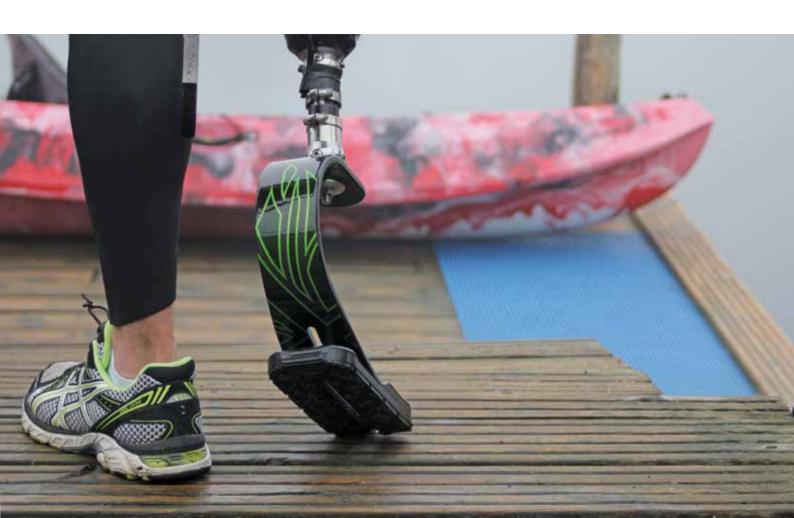

#### Optimal Energy Response and Acceleration

The toe spring provides shock absorption and energy return at toe-off. From the initial ground strike, the C-shaped curve compresses and stores energy, ready to release and provide forward progression into the user's next step.

#### Stability and Safety

With a unique toe and heel design, users feel balanced and secure when standing, running or making a change of direction. An integral traction sole helps to provide safety and sure-footedness for all activities.

#### **Protection and Compliance**

This unique design also ensures shock absorption for comfort and protection, making this the ideal solution for an active user. The split toe improves ground compliance, provides flexibility and enhances comfort on uneven ground. The BladeXT design reinforces natural running motion, reducing strain on the back and shoulders.

#### Lightweight and Strong

Lightweight and strong, BladeXT's design won't slow you down. Made from e-carbon, the BladeXT is durable, yet flexible and responsive.

#### Waterproof

Been somewhere muddy? Simply rinse the BladeXT with fresh water after use in abrasive environments and after exposure to salt or chlorinated water.



#### **Features**

- Twin toe and heel springs
- Split toe design
- Traction sole
- Heel stiffness wedge
- Waterproof
- Field replaceable heel spring kit







#### **Specifications**


Max. User Weight: 166kg (365 lb)

Activity Level: 4
Spring range: 1-9

Component Weight: 850g (1 lb 14 oz) cat. 6
Build Height: 250mm (female pyramid)
230mm (male pyramid)

 $(9^{1}/_{16} - 9^{27}/_{32}")$ 

### **Build Height**

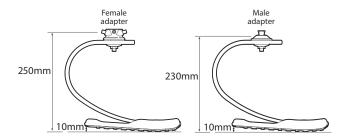



Suitable for submersion

#### Order Example

Product Code	Spring Set	Sole	Adapter Type
BXTAG	3	S	F

Foot example: BladeXT, spring rating 3, sport sole, female adapter

Product Code	Spring Set	Sole	Adapter Type
BXTAG	4	S	М

Foot example: BladeXT, spring rating 4, sport sole, male adapter

#### **Adaptors**

Proximal Adaptor Male 189427
Proximal Adaptor Female 189527

Replacement Heel Spring Kits:

 1, 2 (soft)
 539083SS

 3, 4, 5 (medium)
 539083MS

 6, 7 (firm)
 539083FS

															L
Jogger Occasional to frequent recreational running			اج	Jser We	eight									7	Į
recreational fullilling	<u>r</u>	:≥	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166	kg	ı	
Runner Frequent distance running and endurance training	User	Act	100-115	116-130	131-150	151-170	171-195	196-220	221-255	256-285	286-325	326-365	lbs		
	Jogger	4			3	4	5	6		8	9	9	et e	l	
Important: For higher impact users, do not exceed the weight limit for individual springs.	ooggo.				· ·					Ŭ	· ·	Ŭ	_ 5 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0	ı	
	Runner	4	2	3	4	5	6	7	8	9	9	=	Foot spring set		
	· opgo.														ı



An annual visual inspection is recommended. Check for visual defects that may affect proper function. Maintenance must be carried out by competent personnel. Before carrying out any new activities of daily living, please check with your clinician whether specific training is required.

800 548 3534 | info@blatchfordus.com Blatchford Inc., 1031 Byers Road, Miamisburg, Ohio 45342, USA.

@blatchfordUS | blatchfordus.com

